

Anti Human Organic Anion Transporter 3 (OAT3) Polyclonal Antibody


To eliminate the drug, xenobiotics, a variety of endogenous substances, and their metabolites out of the body, specific membrane proteins named transporters are required. There are two major pathways for the elimination, one of which is hepatic one through bile, and another is renal one to urine. The transporter fall into various transport systems by the transportative substrate. In particular, organic ion transporter family is comprised of organic anion transport family (OAT), organic cation transport family (OCT), OCTN/carnitine transport family, and OAT are multispecific organic anion transporters, the substrates of which include a lot of both endogenous and exogenous anions.

Human organic anion transporter 3 (OAT3) encodes a 543 amino acid residue protein, and which is predicted 12 putative membran-spanning protein. OAT3 is expressed in the kidney, brain, and skeletal muscle. OAT3 mediated the uptake of organic anions, such as PAH (p-aminohippurate), ochratoxin A and estrone sulfate, prostaglandin E₂ (PGE₂), and cimetidine.

This has been proved to be useful for immunohistochemistry .

Package Size	25 µg (250 µL / vial)
Format	Rabbit polyclonal antibody 0.1 mg/mL
Buffer	Block Ace as a stabilizer, containing 0.1%Proclin as bacteriostat
Storage	Store below -20 Once thawed, store at 4 °C. Repeated freeze-thaw cycles should be avoided.
Purification method	This antibody was purified from rabbit serum immunized with synthesized C-end peptide of human OAT3 by peptide affinity chromatography.
Working dilution for immunohistochemistry	1-5 µg/mL
HGNC Name	SLC22A8 (solute carrier family 22A8)

*HGNC: Human Gene Nomenclature Committee

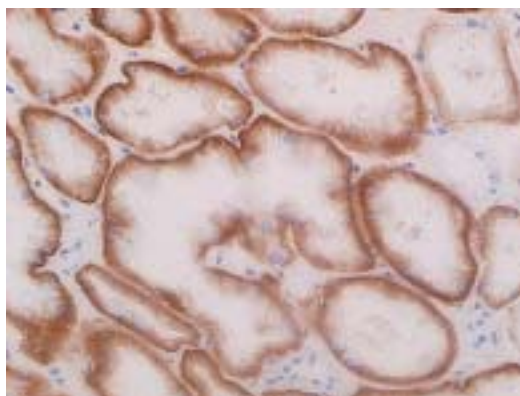

Human Kidney (frozen section)

Anti Human Organic Anion Transporter 3 (OAT3) Polyclonal Antibody

【Reference】

- 1 Kusahara H., Sekine T., Utsunomiya-Tate N., Tsuda M., Kojima R., Cha S.H., Sugiyama., Kanai Y.and Endou H.:Molecular cloning and characterization of new multispecific organic anion transporter from rat brain. *J.Biol.Chem.*274 (19) 13675-13680, 1999
- 2 Sekine T., Cha S.H., Kanai Y.and Endou H.:Molecular biology of multispecific organic anion transporter family (OAT family). *Clin.Exp.Nephrol.*3.237-243,1999
- 3 Sekine T., Cha S.H.and Endou H.:The multispecific organic anion transporter (OAT) family. *pflugers Arch-Eur.J.Physiol.*440.337-350,2000
- 4 Cha S.H., Sekine T., Fukushima J., Kanai Y., Kobayashi Y., Goya T.and Endou H. :Identification and Characterization of Human Organic Anion Transporter 3 Expressing Predominantly in the Kidney. *Mol Pharmacol.* 59.1277-1286,2001

Manufacturer



7-1-14 Minatojimaminami-machi, Chuo-ku, Kobe, Japan 650-0047

Telephone: +81-78-306-0295 FAX:+81-78-306-0296

URL:<http://www.transgenic.co.jp> techstaff@transgenic.co.jp